


- The product 0803922 is a modern single-channel electrocardiograph designed for reliable and easy recording of electrocardiographic signals. It features a thermal ECG waveform printing system, smart warnings for Lead off (electrode disconnection) and Lack of paper, a menu in Greek and English, ECG recording in AUTO (automatic) and Manual (manual) modes, and an ergonomic design and simple operation, suitable for healthcare professionals.

Electrocardiograph 1 Channel with Printer



Ref. 0803922

2 years
guarantee

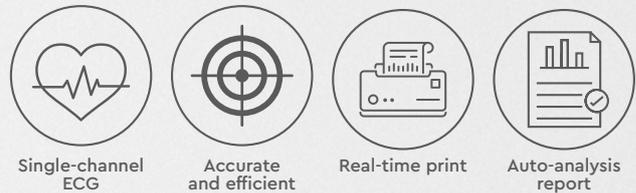

Features:

- Uses advanced digital signal processing technology, offering high-quality ECG waveforms through power, baseline, and EMG filters.
- It incorporates a high-resolution 8 dot/mm thermal printing system, with no need for adjustments, with a recording frequency of up to 150Hz.
- Provides clear and accurate HCG waveforms in real time and continuously, with commentary (electrode points, sensitivity, paper speed, filter status, etc.).
- In AUTO mode, recording is completed with a single touch, significantly increasing efficiency.
- MENU available in Greek and English
- With optimal DC usage, it offers up to 7 hours of standby time, continuous printing for up to 4 hours, and the ability to record up to 150 ECG waveforms.

Package included:

- Lead wire
- Set of limb clips
- Set of chest electrodes
- Printed paper
- Power cord
- Ground wire
- Instruction manual

- Product Dimensions (LxWxH): (31.5x21.5x7.7)cm
- Net weight: 1.5Kg